

Security you can rely on.

Multilink Automix

NEW "Easy Clean-Up" Formula

Multilink[®]
Automix
UNIVERSAL ADHESIVE CEMENT




ivoclar
vivadent
passion vision innovation

Strong


Multilink Automix features patented, hydrolytically stable acidic monomers generating high immediate bond strengths and durable adhesion for long-term security. Within 10 minutes, Multilink Automix establishes a reliable and long-lasting bond to a large variety of restorative materials.

Shear Bond Strength Test Results


All-ceramic to Dentin (Immediate)

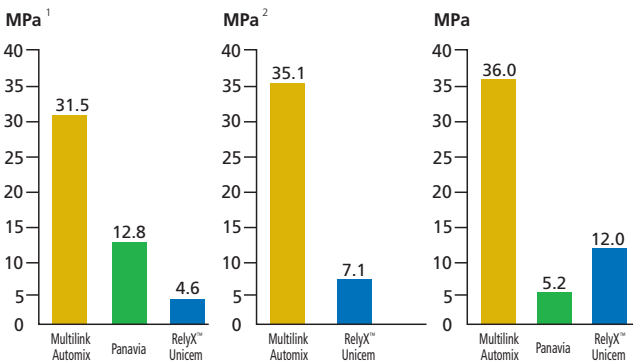

Zirconia to Dentin (Immediate)

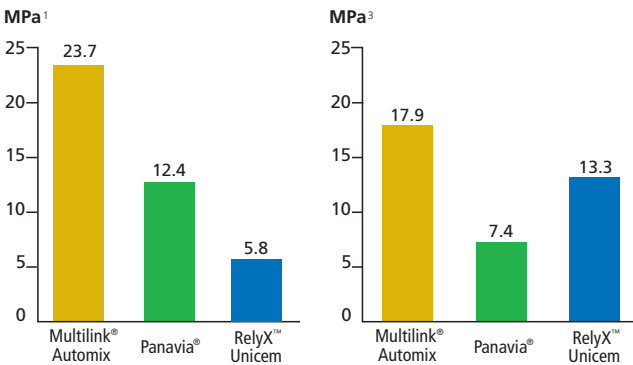

Metal-based to Dentin[†]

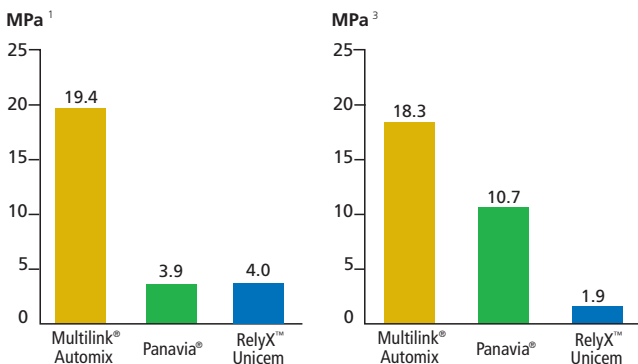

Universal

Multilink Automix is a truly multi-purpose adhesive cement, indicated for cementation of indirect restorations (inlays, onlays, crowns, bridges, and posts) made of all types of restorative materials including:

Metal and Metal-Ceramic Restorations requiring fast and reliable self curing (chemical curing)



IPS InLine®/IPS InLine® PoM

Oxide Ceramics and Fiber-Reinforced Composite Restorations requiring dual-curing initiators due to restricted or hampered access to light (e.g. Zirconium)



IPS e.max® ZirCAD/Ceram



FRC Postec Plus®

All-Ceramic Restorations requiring high-immediate bond strength, good marginal adaptation and effective sealing of the dentin.



IPS Empress®/IPS e.max®

Precious Alloy Restorations requiring long-term reliable results.



Brite Gold® XH

Applied per manufacturer's instructions; Light-cured with high output quartz halogen curing light (>1,000 mW/cm²); Full report available upon request

¹Ivoclar Vivadent Applied Testing Center Amherst, NY, February 2005

²James Simon, DDS, Med Waldemar deRijk, DDS, Ph.D. Director of Esthetic Dentistry and Clinical Research University of Tennessee, February 2005

³Dr. Carlos Munoz Chairman of Restorative Dentistry University at Buffalo, February 2005

[†]10 minute bond strength results; self-cure.

Multilink Automix is a trademark of Ivoclar Vivadent, Panavia is a trademark of Kuraray Co. Ltd and RelyX Unicem is a trademark of 3M ESPE.

Fast, Easy Application

Multilink Automix is a universal adhesive resin cement designed for fast, easy application.


One-Step Primer



A quick 15 second application of the self-etching, self-curing Multilink Primer seals the dentin, provides good marginal adaptation and high immediate bond strength.

Automix Delivery

- Exact dosing with less waste (Dispense only the amount of material required)
- Consistent mix every time
- No hand mixing
- No activation of capsules
- No additional application tools or devices

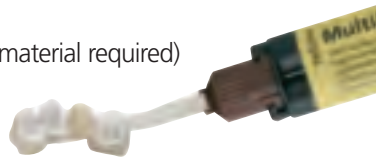


NEW Easy Clean-up Formulation

Multilink Automix is now available in a new "Easy Clean-up" formula for situations where pre-curing excess cement is preferred.



The new "Easy Clean-up" formulation allows excess material to be light-cured for 1-2 seconds per quarter surface (mesio-oral, disto-oral, mesio-buccal, disto-buccal) transforming the material to a gel-like consistency for fast, easy removal.

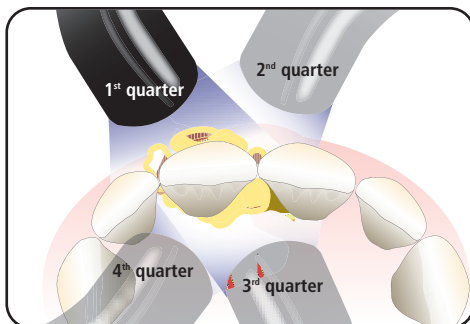

Successful

Multilink Automix is clinically proven by various clinical studies*. The simplified procedure of Multilink Automix in combination with a proven chemistry provides high bond strength values, good marginal adaptation and minimal post-operative sensitivity.



before



prime



seat



clean up



cure



after

*In-vivo studies:

Klinische Langzeituntersuchung von mit Multilink befestigten Porzellanverblend-Restorationen (IPS d.SIGN®) gegenüber IPS Empress®II-Vollkeramik-Restorationen; **Prof G Arnetzel, Universität Graz, Austria (2003)**

Pilot Clinical Trial on Modified Lithium Disilicate Ceramic Crowns; **Prof JA Sorensen, Pacific Dental Institute, OR, USA (2004)**

Clinical Evaluation of an Experimental Dental Ceramic Material for Anterior and Posterior Crowns

Dr Nathanson, Boston University, Boston, MA, USA (2004)

Klinische Evaluation von Multilink an Inlay- und Onlay-Restorationen

Dr A Peschke, Ivoclar Vivadent AG, Schaan, Liechtenstein (2005)

Delivery Forms



System Packs

Contents:

- 1 – Multilink Automix syringe (9g)
- 1 – Primer A (3g)
- 1 – Primer B (3g)
- 1 – Multilink Primer mixing well
- 1 – Monobond Plus (5g)
- 1 – OptraStick test pack
- 15 – Automix mixing tips
- Various accessories

Multilink Automix “Easy Clean-up” System Pack

- 627471 – Transparent
- 627473 – Yellow
- 627472 – White Opaque



Refills

Contents:

- 1 – Multilink Automix syringe (9g)
- 15 – Automix mixing tips

Multilink Automix “Easy Clean-up” Refill

- 615216 – Transparent
- 615217 – Yellow
- 615218 – White Opaque



- 576825 Multilink Primer A/B Refill, (2 x 3g)
- 613626 Multilink Primer A Refill (3g)
- 613627 Multilink Primer B Refill (3g)

626221 Monobond Plus (5g)



592225 Multilink Automix Mixing Tip Refill (15/pk)


603018 Multilink Root Canal Tip Refill (5/pk)

592435 Multilink Applicator Brush Refill (50/pk)

612600 OptraStick Refill (50/pk)

Shades

Multilink Automix is available in 3 shades to meet the clinical indication and esthetic demands.

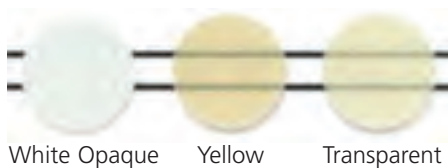

White Opaque Yellow Transparent

Multilink[®]
Automix
 UNIVERSAL ADHESIVE CEMENT

100% CUSTOMER SATISFACTION
 GUARANTEED!

ivoclarvivadent.com

©2010 Ivoclar Vivadent, Inc.
 Multilink, IPS InLine, IPS e.max, FRC Postec Plus, IPS Empress, Brite Gold are registered trademarks of Ivoclar Vivadent, Inc. Clinical photography courtesy of Dr. Robert Ritter, Palm Beach Gardens, Florida
 594992 Rev. 4 12/2010

United States

175 Pineview Drive
 Amherst, NY 14228
 716.691.0010
 800.533.6825
 Fax: 716.691.2285
mail@ivoclarvivadent.com

Canada

2785 Skymark Avenue, Unit 1
 Mississauga, Ontario L4W 4Y3
 905.238.5700
 800.263.8182
 Fax: 905.238.5711
mail@ivoclarvivadent.com

ivoclar
vivadent
 passion vision innovation